
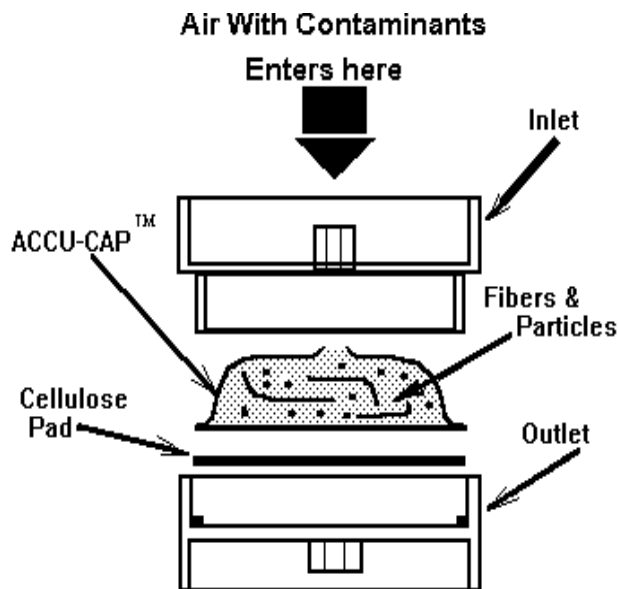


- Prevents loss of sample when removing filter from cassette
- All contaminant is captured in capsule
- Capsule is static dissipative
- Permits in situ extraction of contaminant in capsule

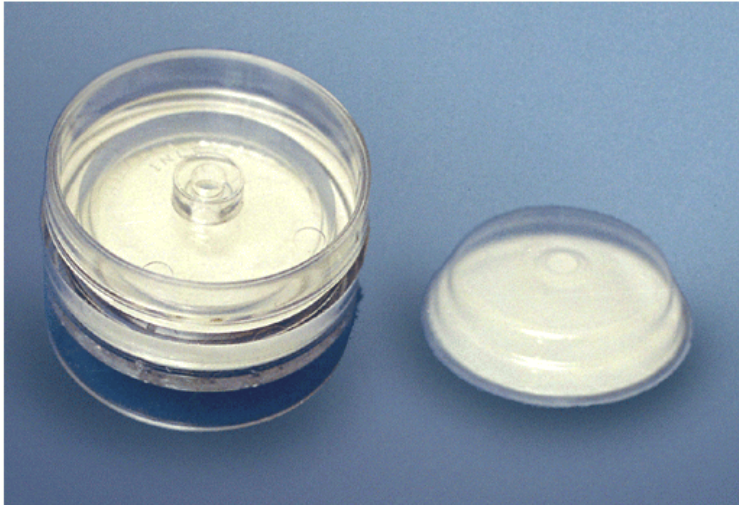

The **ACCU-CAP™** provides for a convenient method for the gravimetric sampling of fibrous dusts. You may have experienced problems when sampling for wood or paper fibers. After sampling, the fibers tend to cling to the inside surfaces of the sampling cassette. Upon opening the cassette to remove the filter for gravimetric analysis, great care must be taken to assure that all the contaminant remains on the filter.

With the use of **ACCU-CAP™**, this difficulty is eliminated; all the fibers are trapped in a removable capsule. The capsule is manufactured from a clear static dissipative plastic, heat sealed to a standard 37mm PVC filter. The capsule fits a 2-piece cassette. Sampling is performed as with a PVC filter; the capsule is then removed and weighed.

Since there is no loss of sample due to electrostatic attraction to the inside surfaces of the sampling cassette, the **ACCU-CAP™** has proven useful in sampling particulate dusts as well.

The **ACCU-CAP™** is available in PVC 5.0µm, 0.8µm; polycarbonate 0.4µm and an all PTFE version.

SPECIFICATIONS

<u>Filter</u>	PVC	Polyvinyl Chloride
	MEC	Mixed Ester of Cellulose
	PC	Polycarbonate
<u>Pore Size</u>	5.0,	0.8, 0.4
<u>Diameter</u>	37mm, 25mm	
<u>Qty/pk</u>	12 or 60	

ORDERING INFORMATION

<u>Description</u>	<u>Qty</u>	<u>Cat#</u>
ACCU-CAP 5.0µm PVC 37mm	12	P-503712-CP
ACCU-CAP 5.0µm PVC 37mm	60	P-503760-CP
ACCU-CAP 0.8µm PVC 37mm	12	P-083712-CP
ACCU-CAP 0.8µm PVC 37mm	60	P-083760-CP
ACCU-CAP 0.8µm MEC 37mm	12	M-083712-CP
ACCU-CAP 0.8µm MEC 37mm	60	M-083760-CP
ACCU-CAP 0.4µm PC 37mm	12	PC043712-CP
ACCU-CAP 0.4µm PC 37mm	60	PC043760-CP

25mm Format

ACCU-CAP 5.0µm PVC 25mm	12	P-502512-CP
ACCU-CAP 5.0µm PVC 25mm	60	P-502560-CP
ACCU-CAP 0.8µm PVC 25mm	12	P-082512-CP
ACCU-CAP 0.8µm PVC 25mm	60	P-082560-CP
ACCU-CAP 0.8µm MEC 25mm	12	M-082512-CP
ACCU-CAP 0.8µm MEC 25mm	60	M-082560-CP
ACCU-CAP 0.4µm PC 25mm	12	PC042512-CP
ACCU-CAP 0.4µm PC 25mm	60	PC042560-CP



4 Kidder Rd. P.O. Box 160 Chelmsford MA 01824-0160

1-800-346-8253 USA and Canada Tel 1-978-256-5450 FAX 1-978-256-8015